
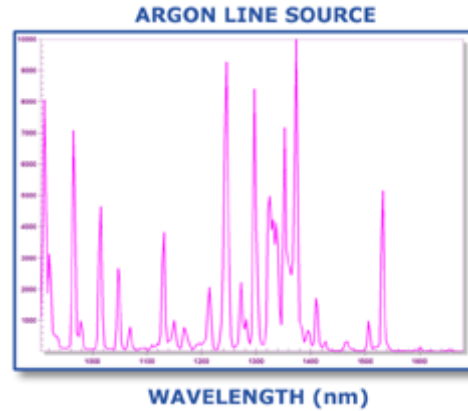


AutoCal

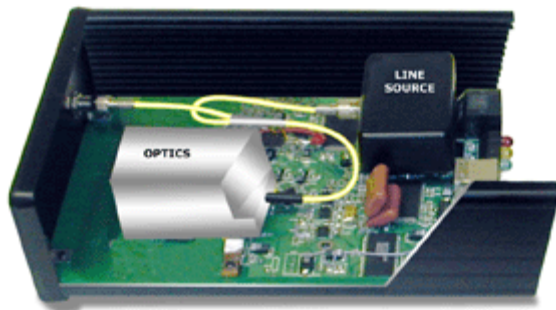

Our AutoCal option provides an automatic wavelength calibration of your spectrometer which guarantees a wavelength accuracy of 0.5nm* or better (typically 0.2*nm or better).

* Dependent on spectral range and slit size

How it works

All of our spectrometers are fiber optic input based. We add a second fiber internally. This second fiber is connected to a line source, Hg/Ar or Xe depending on spectrometer model. The AutoCal feature is totally automated and transparent to the user. The spectrometer recalibrates the wavelength every time a background is saved.

Part Number

AC1 when ordered with spectrometer

AC2 retrofit existing spectrometer

Control Development, Inc., reserves the right to change specifications without notice.

Control Development, Inc.
2633 Foundation Drive, South Bend, IN 46628
Phone: (574) 288-7338 Fax: (574) 288-7339
www.controldevelopment.com
sales@controldevelopment.com

